



Haddii lagugu daweynayo Novolizer®

Novolizer® – way sahlantahay in si sax ah loo isticmaalo

Riix badhanka qaybiyaha si aad u-diyaariso qaadasho mar ah. Marka uu tilmaamuhu cagaar noqdo, xaddiga qaadashada waa ay ku jirtaa oo isku-buufintu waa diyaar. Marka si sax ah la isugu buufiyo wuxu muujiyuhu iska beddelay cagaar oo wuxu noqon guduud. Haddii midabku is-beddeli waayo mar kale isku-buufi.

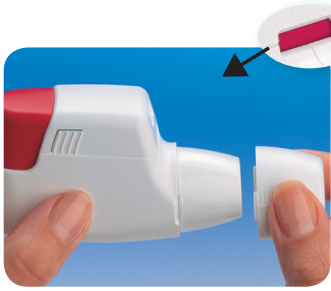

1. Daboolka ka qaad

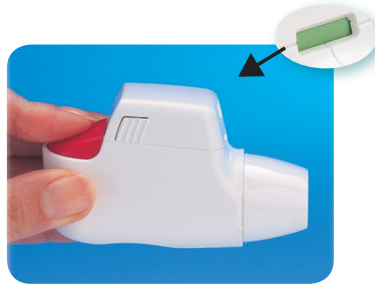


2. Guji badhanka cabbaynta



3. Neefta qaado

Novolizer® – waad garan marka aad si sax ah u samayso

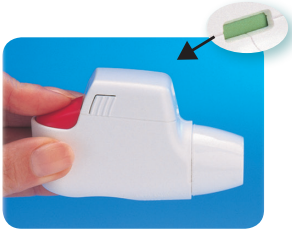

Laba jeer guji

Xaajinta mar-qaadashada ayaa diyaar ah



Cagaar

Qalabka neef-qaadashadu waa u diyaar isticmaalid



Guji

Waxay caddaynaysaa neef-qaadashada saxda ah



Gaduud

Waxay caddaynaysaa neef-qaadashada saxda ah



Dhadhamin

Waxay caddaynaysaa neef-qaadashada saxda ah



Ku iskaan-garee sumadda QR telefonkaaga moobilka si aad u aragto tilmaamaha fiidiyaha ee www.novolizer.se

Novolizer® – xaddi-dawo tiriye muuqanaya

Novolizer® waxay leedahay shaashad muujisa tobanbkii mar-qaadasho ee kasta. Marka tiriyuhu muujiyo 0 waxa la joogaa wakhtigii la bedeli lahaa cajalka.



Booqo websaydyada:


- www.novolizer.se
- www.astmalinjen.se
- www.kolguiden.se

Novolizer®

Qalab neef-qaadasho oo isbedel keenaya

Geli qalabka neef-qaadashada cajal cusub


1. Si sahlan u cadaadi meelaha godgodan ee dhinacyada daboolka, hore u riix daboolka, dabadeedna kor uga soo qaad.

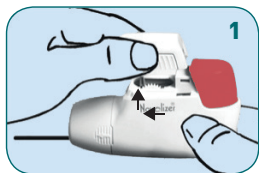


2. Dahaarka almuniyada ah ee ilaalinaya Ka jar galka cajalka oo soo saar cajalka.



3. Cajalada geli isku-buufinta iyada oo shaashadda mar-qaadashadu ujeedo dhinaca hore.

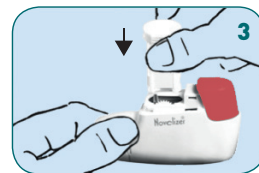

4. Dib ugu celi daboolka adiga oo dhinaca kore ka soo geli-naya oo hoos u cadaadi ilaa aad maqlayso cod yar.

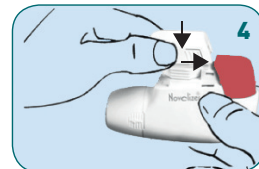

Novolizer® – qalab neef-qaadasho oo isbedel keenaya

Novolizer® waxa loo isticmaalaa daawaynta xiiqda (neefta) iyo jirrooyinka xiritaanka sanbabada “COPD (chronic obstructive pulmonary disease)”.

Novolizer® waxay ku timaadaa taxanayaal:



Novopulmon® Novolizer® (budesonide)
Neef-qaadashada steroid-ka

200 µg, 200 xaddi-dawo
400 µg, 200 xaddi-dawo



Formatris® Novolizer® (formoterol)
Muddo dheer waxtareysa hawo maren ballaariye

6 µg, 60 xaddi-dawo
12 µg, 60 xaddi-dawo



Ventilastin® Novolizer® (salbutamol)
Muddo gaaban waxtareysa hawo maren ballaariye

100 µg, 200 xaddi-dawo

Novolizer® waxay gaaraysaa ilaa 2000 oo neef-qaadasho ama ilaa hal sano. Taas macneheedu waxa weeye in aad adigu iskaa u cabbayn karto qalabka neef-qaadashada, taas oo yaraynaysa qarashka iyo saamaynta cimilada labadaba.

Qorshaha daweyntayda

Novopulmon® Novolizer® – neef-qaadashada steroid-ka

Quwadda: 200 µg

400 µg

Tirada neef-qaadasho/maalintii: Taariikhda bilawga:

Ventilastin® Novolizer® – muddo gaaban waxtareysa hawo maren ballaariye

Quwadda: 100 µg

Tirada neef-qaadasho/maalintii: Taariikhda bilawga:

Formatris® Novolizer® – muddo dheer waxtareysa hawo maren ballaariye

Quwadda: 6 µg

12 µg

Tirada neef-qaadasho/maalintii: Taariikhda bilawga:

Novolizer®

Qalab neef-qaadasho oo isbedel keenaya